


Alpha Transfer System 400

High SWL
Safe lifting

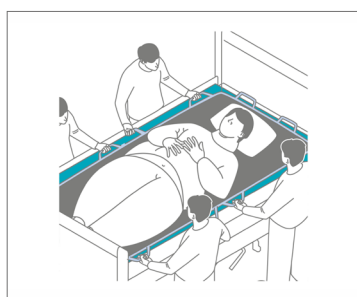

The Alpha Transfer System 400 is a lifting sheet for clients up to 400 kg/881 lb. The person on the sheet can be carried by a maximum of 12 persons. There are four hand-grips on the long sides, two on the head side and two on the foot side. Every handgrip is attached to the opposite grip in the transfer system. The grip covers make holding easier.

Our OnlyU Disposable Sliding Foil can be used on the Transfer System to assist in moving clients on solid surfaces and bases. Friction can be significantly reduced.

Technical specification


Product Name	Art. nr.	Dimensions	Material	Max lifting weight
Transfer System 400	PM-1501	115 x 220 cm / 45,2" x 86,6"	PES with PVC	400 kg / 881 lb



CE (MD) = 400 kg / 881 lb

PETERMANN
part of the Human Care Group

These products are distributed by
Human Care AB for Petermann GmbH.

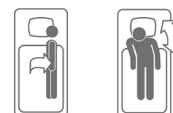

Visit www.humancaregroup.com for the latest product information and contact details to national offices and partners.